


How to de-identify data in FLIP datafile

Raw flip datafile will be like:

xxxxLIVE_DATA_1(or2).TXT

Open the file by Chrome (or Firefox)

Or as below: Drag the file: 'wk_Endoflip_remove_data.html' to Chrome browser

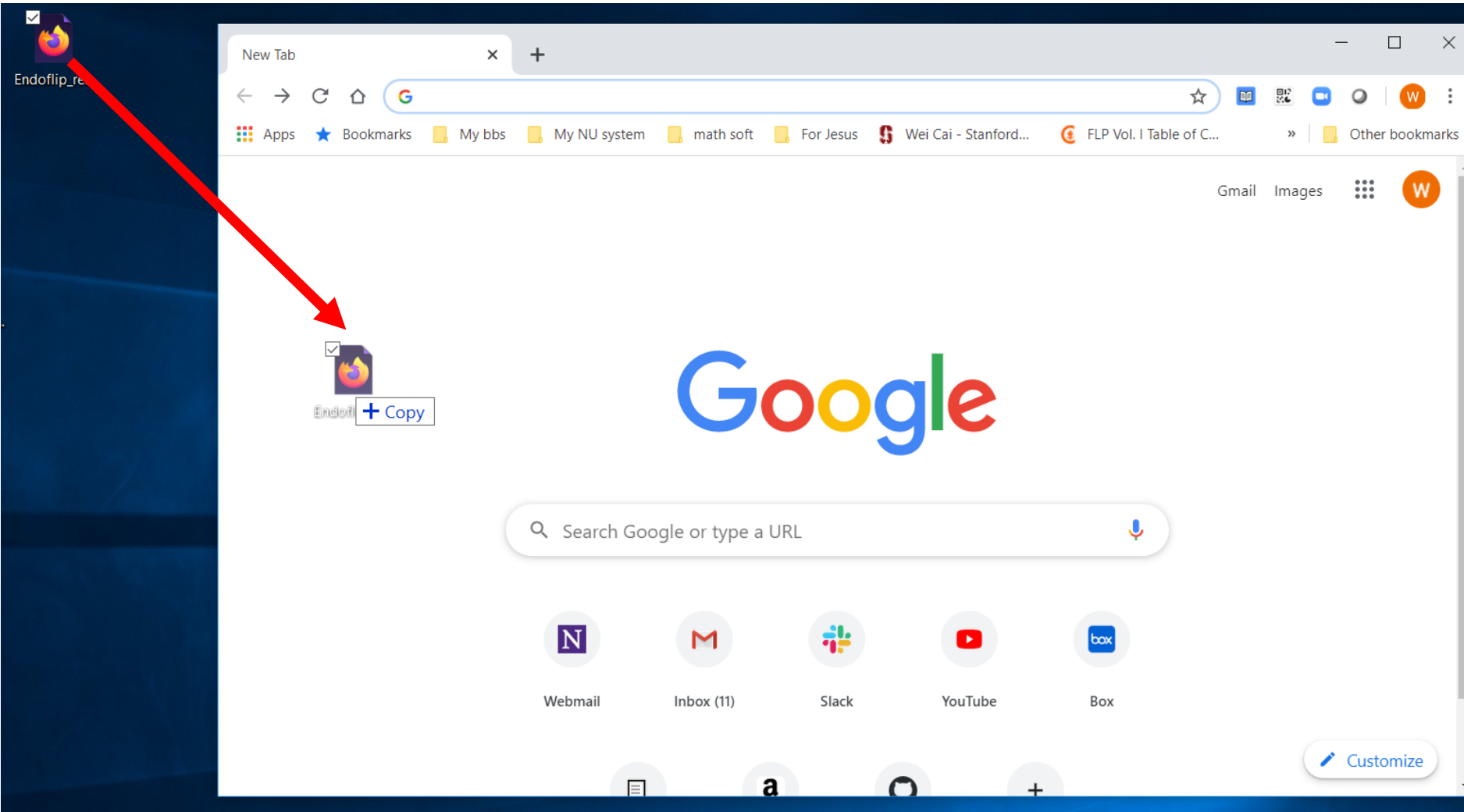

2. Once the file is open, a webpage like below should be shown.
In the webpage: the upper are two text fields; the lower are steps to follow.

Declassify date in Endoflip

File | C:/Users/wkt699/Desktop/Endoflip_remove_date.html

Apps Bookmarks My bbs My NU system math soft For Jesus Wei Cai - Stanford... FLP Vol. I Table of C... Other bookmarks

Declassify: date in Endoflip? Follow the three steps (bottom).

[Click to see HIPAA Policy](#)

Original Flip file (first 5 lines in the file):

Test declassified data (first 5 lines in the file):

Step 1- Select a File to Load:

Choose File No file chosen Load Selected File

Step 2- Declassify the file (i.e. Replacing date):

Original Date: Will be replaced by: Click to test declassified data

Step 3: Save the new declassified file as:

New Filename: Click to save new File

- 1: Click 'Choose File' and you will see window to locate file;
2. Click 'Load selected File', and you will see data in the upper text field.

Declassify date in Endoflip

File | C:/Users/wkt699/Desktop/Endoflip_remove_date.html

Apps Bookmarks My bbs My NU system math soft For Jesus Wei Cai - Stanford... FLP Vol. I Table of C... Other bookmarks

Declassify date in Endoflip? Follow the three steps (bottom).

[Click to see HIPAA Policy](#)

Original Flip file (first 5 lines in the file):

```
Time,Entry No.,Pump
Status,E1DS100*,E2DS100*,E3DS100*,E4DS100*,E5DS100*,E6DS100*,E7DS100*,E8DS100*,E9DS100*,E10DS100*,E11DS100*,E12DS100*,E
13DS100*,E14DS100*,E15DS100*,E16DS100*,BP1,BV,BT1,Catheter Current90,P.Sense

1990:01:01-08:18:13|408,19325,S, 5.38, 5.38, 5.38, 5.38, 5.38, 5.38, 5.38, 5.38, 5.38, 5.38, 5.54, 5.38, 5.38, 5.38,
5.84, 5.38, -0.15, 0, 17.47, 83.52, 0

1990:01:01-08:18:13|508,19326,S, 5.38, 5.38, 5.38, 5.38, 5.38, 5.38, 5.38, 5.38, 5.38, 5.38, 5.54, 5.38, 5.38, 5.38,
5.84, 5.38, -0.15, 0, 17.47, 83.50, 0
```

Test declassified data (first 5 lines in the file):

1

2

Step 1- Select a File to Load:

Choose File test_LIVE_DATA_1.TXT Load Selected File

Step 2- Declassify the file (i.e. Replacing date):

Original Date: 1990:01:01 Will be replaced by: xxxx:xx:xx Click to test declassified data

Step 3: Save the new declassified file as:

New Filename: case1.TXT Click to save new File

- 3: Click 'Click to test declassified data' and you will see lower text field with replaced date
- 4: Click 'Click to save new File', and you will see a new file downloaded. You can rename the file.

Declassify date in Endoflip

File | C:/Users/wkt699/Desktop/Endoflip_remove_date.html

[Click to see HIPAA Policy](#)

Original Flip file (first 5 lines in the file):

```
Time,Entry No.,Pump
Status,E1DS100*,E2DS100*,E3DS100*,E4DS100*,E5DS100*,E6DS100*,E7DS100*,E8DS100*,E9DS100*,E10DS100*,E11DS100*,E12DS100*,E13DS100*,E14DS100*,E15DS100*,E16DS100*,BP1,BV,BT1,Catheter Current90,P.Sense

1990:01:01-08:18:13|408,19325,S, 5.38, 5.38, 5.38, 5.38, 5.38, 5.38, 5.38, 5.38, 5.38, 5.38, 5.54, 5.38, 5.38, 5.38, 5.84, 5.38, -0.15, 0, 17.47, 83.52, 0

1990:01:01-08:18:13|508,19326,S, 5.38, 5.38, 5.38, 5.38, 5.38, 5.38, 5.38, 5.38, 5.38, 5.38, 5.54, 5.38, 5.38, 5.38, 5.84, 5.38, -0.15, 0, 17.47, 83.50, 0
```

Test declassified data (first 5 lines in the file):

```
Time,Entry No.,Pump
Status,E1DS100*,E2DS100*,E3DS100*,E4DS100*,E5DS100*,E6DS100*,E7DS100*,E8DS100*,E9DS100*,E10DS100*,E11DS100*,E12DS100*,E13DS100*,E14DS100*,E15DS100*,E16DS100*,BP1,BV,BT1,Catheter Current90,P.Sense

xxxx:xx:xx-08:18:13|408,19325,S, 5.38, 5.38, 5.38, 5.38, 5.38, 5.38, 5.38, 5.38, 5.38, 5.38, 5.38, 5.54, 5.38, 5.38, 5.38, 5.84, 5.38, -0.15, 0, 17.47, 83.52, 0

xxxx:xx:xx-08:18:13|508,19326,S, 5.38, 5.38, 5.38, 5.38, 5.38, 5.38, 5.38, 5.38, 5.38, 5.38, 5.38, 5.54, 5.38, 5.38, 5.38, 5.84, 5.38, -0.15, 0, 17.47, 83.50, 0
```

Step 1- Select a File to Load:

Choose File test_LIVE_DATA_1.TXT Load Selected File

Step 2- Declassify the file (i.e. Replacing date):

Original Date: 1990:01:01 Will be replaced by: xxxx:xx:xx Click to test declassified data

Step 3: Save the new declassified file as:

New Filename: case1.TXT Click to save new File

case1 (1).TXT


File downloaded shortly

3

4

Another case?

Just click Refresh, and start over...

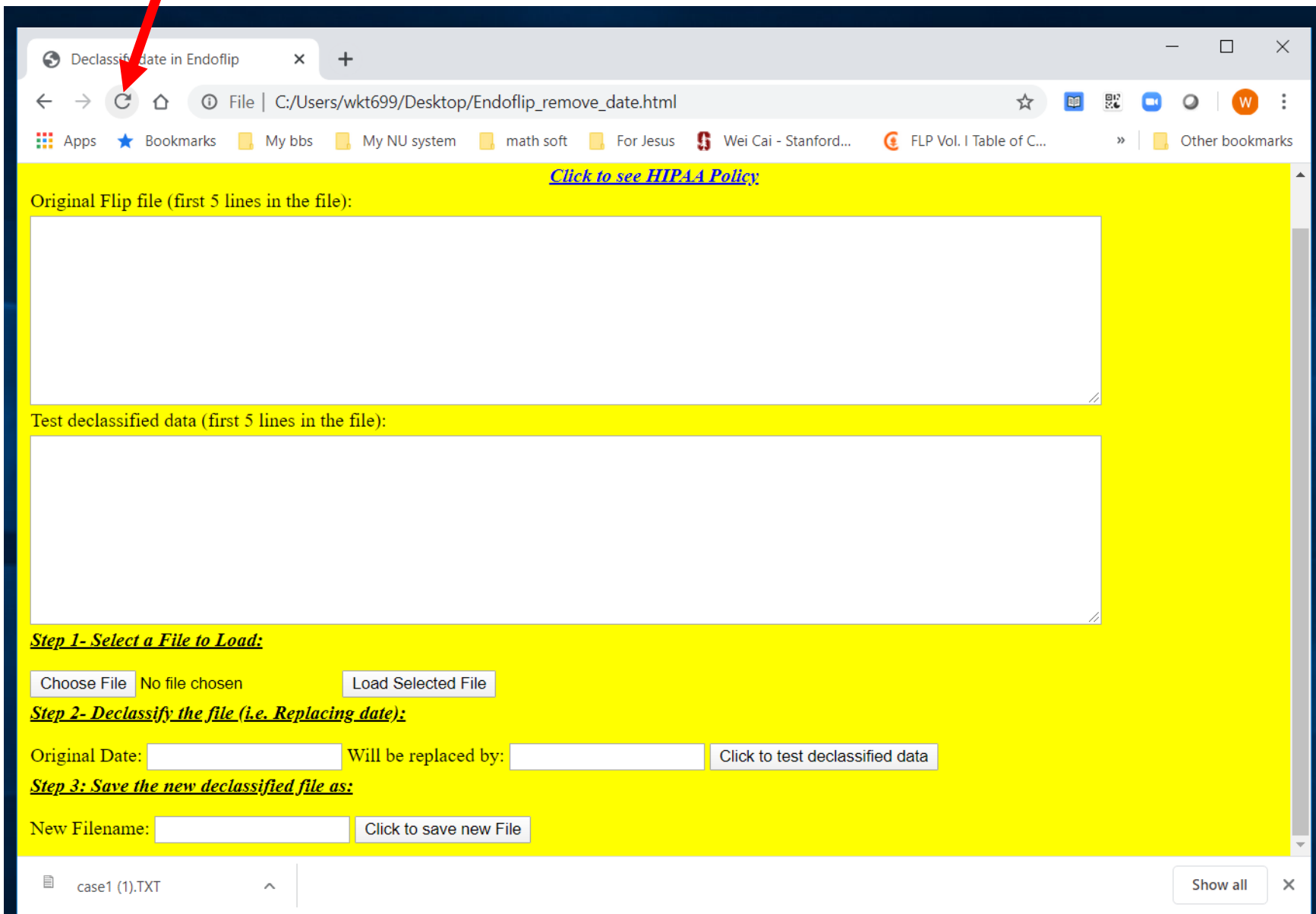

Declassify date in Endoflip

File | C:/Users/wkt699/Desktop/Endoflip_remove_date.html

Apps Bookmarks My bbs My NU system math soft For Jesus Wei Cai - Stanford... FLP Vol. I Table of C... Other bookmarks

[Click to see HIPAA Policy](#)

Original Flip file (first 5 lines in the file):

Test declassified data (first 5 lines in the file):

Step 1- Select a File to Load:

Choose File No file chosen Load Selected File

Step 2- Declassify the file (i.e. Replacing date):

Original Date: Will be replaced by: Click to test declassified data

Step 3: Save the new declassified file as:

New Filename: Click to save new File

case1 (1).TXT Show all